
Validate: A Workflow for Evaluating Performance for GWAS/QTL Tools Using Known-Truth Datasets

*Ann Stapleton, Kurt Michels, and Dustin Landers
University of North Carolina Wilmington*


Abstract

Understanding the effectiveness of Genome-Wide Association (GWAS) and Quantitative Trait Loci (QTL) analytical tools under various situations is crucial to deciding which tools are best given a particular problem. Validate provides a way to return classification and regression performance measures for large amounts of tool outputs generated using known-truth simulations. We also provide solutions for aggregating hundreds or thousands of outputs in to a single folder on the iPlant data store, so that Validate can be used.

Contents

1	Introduction	3
1.1	About the Developers	3
1.2	Why Validate?	3
1.3	Getting iPlant Credentials	4
1.4	Getting Your Application on an API	4
1.5	Where Can I Get These Softwares?	4
1.6	How to Use This Manual	4
2	Running Simulations	6
2.1	How To Run Simulations with rPlant	6
2.2	Using rPlant to Sample From Your Known-Truth Data Sets	9
3	Aggregating with Aggregate	10
3.1	How to Use the GUI Aggregate Tool	10
3.2	How to Use the Aggregate Function in rPlant	16
3.3	How to Use Aggregate on Atmosphere	16
4	Validating with Validate	16
4.1	How to Use Validate on the DE	16

5	Visualizing the Validation	19
5.1	How to Use the Demonstrate R Package	19
5.2	How to Use Demonstrate on the DE	19
5.3	How to Use Demonstrate on Atmosphere	19
6	Acknowledgements	19
7	How to Interpret Performance Measures	19
7.1	AUC	19
7.2	H	19
7.3	KS	19
7.4	TPR	19
7.5	FPR	19
7.6	Accuracy	19
7.7	MAE	19
7.8	RMSE	19

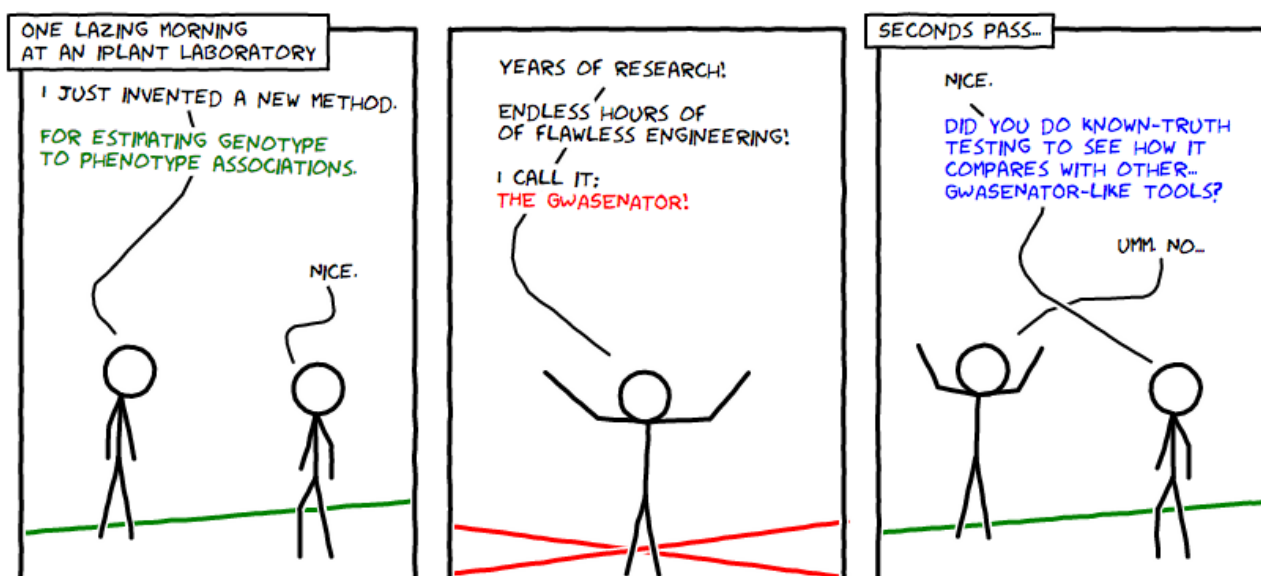

1 Introduction

1.1 About the Developers

Project Lead and Principal Investigator

Dr. Ann Stapleton works at the interfaces, such as the junctions between research and teaching, individual research projects and large collaborative projects, the organization of international meetings and high-school teacher inquiry labs. Most recently, she has worked at the interface between plant biology and software engineering, leading the way to broad methods applicable to evaluating genotype-to-phenotype analytical methods.

Statistical Analyst and Developer

Dustin Landers received a BS from Appalachian State where he learned a lot about survey research and using statistics to solve problems and answer questions. He continued his education at UNCW by studying Applied Statistics. Moving from world of polling and survey statistics to the that of Big Data, he has an interest in bridging the gap between statistics and software engineering and seeing the combined discipline brought to bear on problems once considered impossible.

1.2 Why Validate?

To test known-truth data sets, we originally weren't sure of the best way to approach the problem. We had the idea of using ROC plots. We wanted to look at true positive rates, false positive rates and accuracy. Of course, this was never really feasible because any one simulation run could be atypical, and all of these methods could only analyze a single simulation run. We asked ourselves how could we run lots of simulations through a tool, and test the outputs in a way that gave us insight in to the performance of these GWAS and QTL tools?

We decided that in order to allow lots of simulations to be tested, we needed a method that that allowed iterations over hundreds or even thousands of simulation-runs. This birthed the idea of Validate, which is a tool that returns these sorts of measures of tool success for massive amounts of outputs.

But how do we run all these simulations? How do we store them all in a single location so we can calculate these measures? How do we visualize them?

These are the problems we sought to solve. For each problem, we have our own solution. We also left them divided. This way you can decide whether to use our solution or to use your own if you have a unique scenario.

1.3 Getting iPlant Credentials

1.4 Getting Your Application on an API

1.5 Where Can I Get These Softwares?

1.6 How to Use This Manual

Evaluating a GWAS/QTL tool with *Validate* is essentially four major steps after your application is installed on either the iPlant Foundation API or the Agave API:

→ Section 2

1. **Running the Tool With Simulations As Inputs.** This step involves deciding what simulations to use and then iterating over those simulations and submitting them as job requests through the API. We recommend using some sort of scripting method that you are comfortable with. At this point, we don't have a standalone application for submitting jobs. We use *rPlant*, which is freely available R package that allows you to connect with iPlant's API layer to submit job requests.

→ Section 3

2. **Aggregating the Outputs into a Single Folder.** This part is a bit of a logistical exercise. You need to put all the tool outputs (from Step 1) that you want to be analyzed using the same known-truth metadata in to aggregate folders. For example, say we are using simulations that are generated using varying levels of heritability. Since the varying heritability values in essence produce different SNP effects, we need to essentially run *Validate* three separate times. *Validate* requires an input on an entire folder, and then iterates over those tool outputs. So the first step here, is to decide how many different runs we need to do and then create that many folders. We provide a GUI tool, *Aggregate* that allows you to select files from multiple folders on your iPlant data store (multiple runs of a single tool) and move those (or aggregate them) in to a single folder.

→ Section 4

3. **Running *Validate* on the Aggregated Folder.** Once you have all your outputs in aggregated folders. You can simply log in to the iPlant Discovery Environment and run *Validate* on that folder. You must have at least two columns with header names in the outputs: The name of the SNP column, and the name of threshold column (such as P-value). Further, if you wish the get back effect size estimation errors, you must also know the column on the estimated effect size column (for example, PLINK's is BETA). Text files with known SNPs and effect sizes must also be included. Once you submit *Validate*, you will receive a notification when it is completed. The *Validate* output will be columns of performance measures for each tool output.

→ Section 5

4. **Running *Demonstrate* on the *Validate* Output.** Finally, once you have outputs from *Validate* (which may be more than one), you


need to combine them and test the differences in your measures between your simulation's parameters. Demonstrate is an R package (to be installed on the iPlant DE) that quickly combines all your results files in a folder you specify and returns sciplot factorial graphs for parameters you specify. It also returns the combined results file so that you can perform your own analyses.

2 Running Simulations

2.1 How To Run Simulations with rPlant

Step 1-A) Install R and rPlant

We will first show you how to batch run simulations with rPlant. First, we are using the latest R version, 3.0.2. When you open the R console, it should look something like this (we are using a Mac, so yours may look slightly different, but it will function essentially the same).

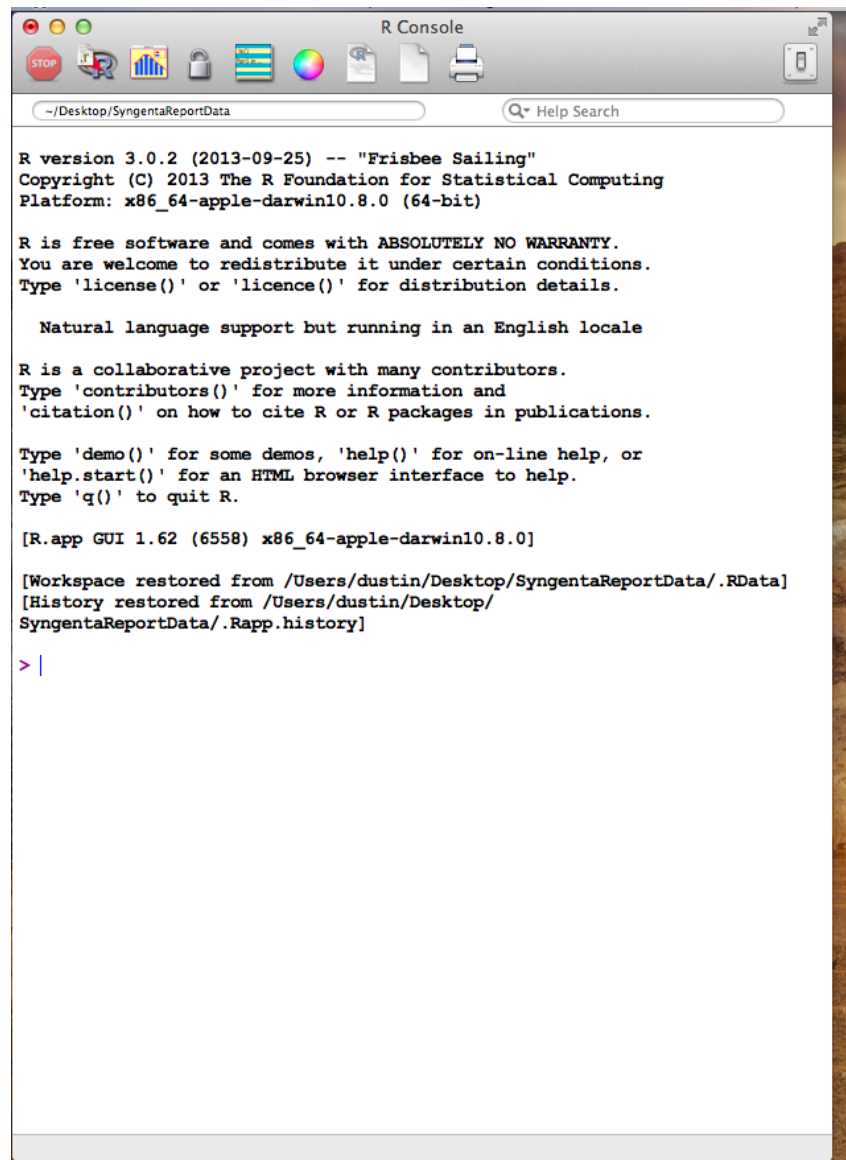

```
R version 3.0.2 (2013-09-25) -- "Frisbee Sailing"
Copyright (C) 2013 The R Foundation for Statistical Computing
Platform: x86_64-apple-darwin10.8.0 (64-bit)

R is free software and comes with ABSOLUTELY NO WARRANTY.
You are welcome to redistribute it under certain conditions.
Type 'license()' or 'licence()' for distribution details.

Natural language support but running in an English locale

R is a collaborative project with many contributors.
Type 'contributors()' for more information and
'citation()' on how to cite R or R packages in publications.

Type 'demo()' for some demos, 'help()' for on-line help, or
'help.start()' for an HTML browser interface to help.
Type 'q()' to quit R.

[R.app GUI 1.62 (6558) x86_64-apple-darwin10.8.0]

[Workspace restored from /Users/dustin/Desktop/SyngentaReportData/.RData]
[History restored from /Users/dustin/Desktop/
SyngentaReportData/.Rapp.history]

> |
```

If you haven't already installed rPlant, we recommend that you do that. After that, you need to validate your iPlant credentials so that you can access your simulation files. Then, you will need to get a list of the names of your simulations files and iterate over them in a for loop, submitting jobs to your installed application.

Step 1-B) Consider the strategy you want to use. Do you want to use all your known-truth data sets, or just a sample? And if a sample, should you take a stratified sample?

→ These are mainly just example codes. Although you are welcome to copy ours.


The following code is to just an example of how the submission may go assuming your known-truth data sets are in something like a raw format, where each known-truth data set is stored in a single file. rPlant assumes you already include your iPlant data store username in the files path, so that "iPlant/home/username/simulations" is just written as "simulations" or "iPlant/home/username/analyses" is just written as "analyses". This is, of course, because validating your credentials allows you to not have to keep typing you user information over and over again.

```
> install.packages("rPlant")
> require(rPlant)
> Validate("username", "password")
> mydir <- ListDir("simulations")
> for (i in 1:nrow(mydir)) {
+   SubmitJob(job.name=NULL,
+             application="myGWASApplication",
+             file.path="simulations", file.list=list(mydir[i,1]))
+   print(i)
+ }
```

Step 1-C) Find out the name of your application on the API, and the inputs required.

You need to make sure of the specific application name assigned to your tool in the API and of any special parameters or inputs required in order to run. For example, PLINK requires a PEDMAP format. The PEDMAP format is essentially two files, the .PED and the .MAP file for each known-truth data set.

For example, our simulation files are stored in a shared folder, where each odd file number is a .PED and each even is a .MAP. Each odd is a .PED file and the subsequent file number is a .MAP like so:

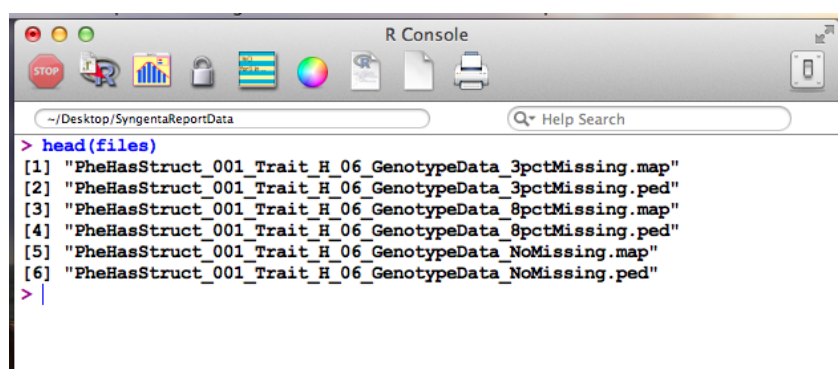

As can be seen, for our particular structure we would need to submit two files for each run. So let's walk through it. We first can list all the apps in order to find ours (which we have already installed), so that we can get the exact name that we need to use for the SubmitJob function.

```
> ListApps()
```

```

[33] "GenomeSelection-0.0.1u2"
[34] "GMAP_stampede-121212u1"
[35] "GSNAP_lonestar-121212u1"
[36] "GSNAP_stampede-121212u2"
[37] "head-stampede-5.97u2"
[38] "idbaUD-1.0.0u2"
[39] "interproscan-5.44.0u1"
[40] "macs-ranger-1.4-1.4u1"
[41] "mafft-7.113u1"
[42] "mafft-lonestar-6.864u1"
[43] "mafftDispatcher-1.0.13100u1"
[44] "MergeG2P-0.0.1u1"
[45] "metagenemark-1.00u3"
[46] "metaphlan-lonestar-1.6.0u4"
[47] "MLMM-0.0.1u1"
[48] "MrBayesmpi_basic-3.2.1u1"
[49] "Muscle-3.8.32u4"
[50] "muscle-lonestar-3.8.31u2"
[51] "newbler-2.6.0u1"
[52] "NPUTE-0.0.2u1"
[53] "NumericalTransform-0.0.1u1"
[54] "oases-0.2.08u1"
[55] "philip-dna-parsimony-lonestar-3.69u2"
[56] "philip-protein-parsimony-lonestar-3.69u2"
[57] "plink-1.07u1"
[58] "prodigal-1.0.0u1"
[59] "quicktree-dm-lonestar-1.1u2"
[60] "quicktree-tree-lonestar-1.1u2"
[61] "raxml-lonestar-7.2.8u1"
[62] "ray-2.2.0u1"
[63] "scarf-1.00u1"
[64] "soapdenovo_trans-1.0u1"
[65] "soapdenovo-1.05u1"
[66] "soapdenovo-2.04u1"
[67] "STRUCTURE-2.3.4u2"
[68] "STRUCTURE-2.3.5u1"
[69] "STRUCTURE2TASSEL-0.0.1u1"
[70] "TASSEL (GLM) -0.0.1u1"
[71] "TASSEL (MLM) -0.0.1u1"
[72] "TASSEL4-GLM-0.0.1u1"
[73] "tasseldDispatcher-1.0.13350u1"
[74] "TNRS4GWAS-0.0.2u1"
[75] "trinity_lonestar4-20130814u1"
[76] "trinity_normalize_by_kmer_coverage_lonestar4-20130814u1"
[77] "velvetg-1.2.07u2"
[78] "velveth-1.2.07u1"
[79] "XYPLOT-0.0.1u1"
>

```

We can see that the exact name of PLINK on the Foundation API is "plink-1.07u1". We will revisit the earlier explanation to see how we would submit a PLINK job using the PEDMAP format.

Step 1-D) Use the short R script to select known-truth data sets and submit a job run for each one using your application.

```

> install.packages("rPlant")
> require(rPlant)
> Validate("username", "password")
> mydir <- ListDir("simulations")
> for (i in 1:nrow(mydir)) {
+   SubmitJob(job.name=NULL,
+             application="plink-1.07u1",
+             file.path="simulations", file.list=list(mydir[i,1],
+             mydir[i+1,1]))
+   print(i)
+ }

```

This should iterate over each simulation in our data set and submit a PLINK job to the iPlant hardware. Since iPlant uses a queuing system,

this effectively distributes your workload over many computer nodes as needed. This is obviously helpful in achieving a high volume of runs, even outside of known-truth testing. Now you know how!

2.2 Using rPlant to Sample From Your Known-Truth Data Sets

Let's say, for example, that we had 600 total known-truth data sets: Our known-truth datasets are generated from hundred different genotypes, with three different types of simulated missing values such as 0%, 3%, and 8%. In addition, we have data generated with and without population structure (meaning the stochastic processes generated the data have either a single population mean or varying group means).

Say that we didn't want to run all of them, but maybe just half of them. In many cases, you may have substantially more than 600 and thus sampling makes much more sense then it does in this example.

We can use some scripting in R and rPlant to accomplish the same thing we did earlier but with some random sampling. Here is the code we used to sample 300 known-truth PEDMAPs for PLINK.

```
> set.seed()
> require(rPlant)
> Validate("username", "password")
> mydir <- ListDir("simulations")
> genos <- matrix(nrow=1200, ncol=1)
> break.data <- function(x, y=2) {
+   broken <- unlist(strsplit(x, '-', TRUE))
+   return(broken[y])
+ }
> mydir <- mydir[, 1]
> for (i in 1:length(mydir)) {
+   genos[i] <- break.data(files[i], 2)
+   gens <- na.omit(unique(genos))
+ }
> oddsamp <- seq(1, 11, 2)
> totest <- list()
> for (i in 1:length(genoa)) {
+   myfiles <- list()
+   for (j in 1:length(files)) {
+     if (breakdat(files[j], 2) == genos[i]) {
+       myfiles <- append(myfiles, files[j])
+     }
+   }
+   myfiles <- myfiles[sample(oddsamp, 2, FALSE)]
+   myfiles <- unlist(myfiles)
+   for (k in 1:length(files)) {
+     for (p in 1:length(myfiles)) {
+       if (files[k] == myfiles[p]) {
+         totest <- append(totest, k)
+       }
+     }
+   }
+ }
> totest <- unlist(totest)
> for (i in 1:length(totest)) {
+   SubmitJob(job.name=NULL,
+             application="plink -1.07ul", file.list=
+               list(files[totest[i]], files[totest[i]+1]),
+               file.path="simulations")
+   print(i)
+ }
```

Every situation is unique in this case, which is part of the reason we have yet to develop a graphical user interface method for submitting jobs. All

that we did in the code above, is demonstrate how rPlant may be used to break up some files that had naming conventions that were separated by underscores. If you look back to the earlier image that showed what our known-truth data sets looked like, you would see that phenotype, heritability information, and population structure are all in the naming format for each data set.

This is an important piece of information, because without an appropriate naming structure, there would be no way for us to keep metadata on the known-truth data sets throughout this testing process.

We break up the names to grab the pieces about phenotype information in order to essentially stratify-sample our data sets by genotype. We then have six data sets for any one genotype (three missing values types and two population structure types), and then we just pick two at random.

This insures that we include all of our genotypes in our runs, but there are simpler ways to do this, that basically just random sample from all of our data sets. The following code again assumes you are using PEDMAPs to run PLINK, but that you just want to sample a random 300 (as opposed to a stratified-random sample).

```
> set.seed()
> require(rPlant)
> Validate("username", "password")
> mydir <- ListDir("simulations")
> mydir <- mydir[,1]
> oddsamp <- seq(1,599,2)
> totest <- sample(oddsamp, 300, FALSE)
> for (i in 1:length(totest)) {
+   SubmitJob(job.name=NULL,
+             application="plink-1.07ul", file.list=
+             list(mydir[totest[i]], mydir[totest[i] + 1],
+             file.path="simulation")
+ }
```


3 Aggregating with Aggregate

3.1 How to Use the GUI Aggregate Tool

The second major step is accomplished once all the runs have finished. In our example, depending on the resources available, PLINK can finish in a range from ten seconds to two and half minutes per run (not including the queuing process). This is non-linear as well, because many will be running simultaneously on separate nodes while some remain in queues. Its hard to gauge exactly how long your tool will take, but it will almost definitely run faster on multiple nodes this way.

Step 2-A) Log in to the DE and make sure all runs are complete.

Let's say that we have already run all of our runs. Then we would want to log in to the iPlant Discovery Environment to make sure all the folders are in our analyses folder. As can be seen, each run of PLINK with each distinct known-truth data set has been put in to each folder. If we were to open up each folder, we would see a set of files, of which only one we will want to put through Validate. In the case of this PLINK run, you need to put each file ending in .qassoc (the standard PLINK output filetype) in to a single folder, so that we can run Validate on that folder.

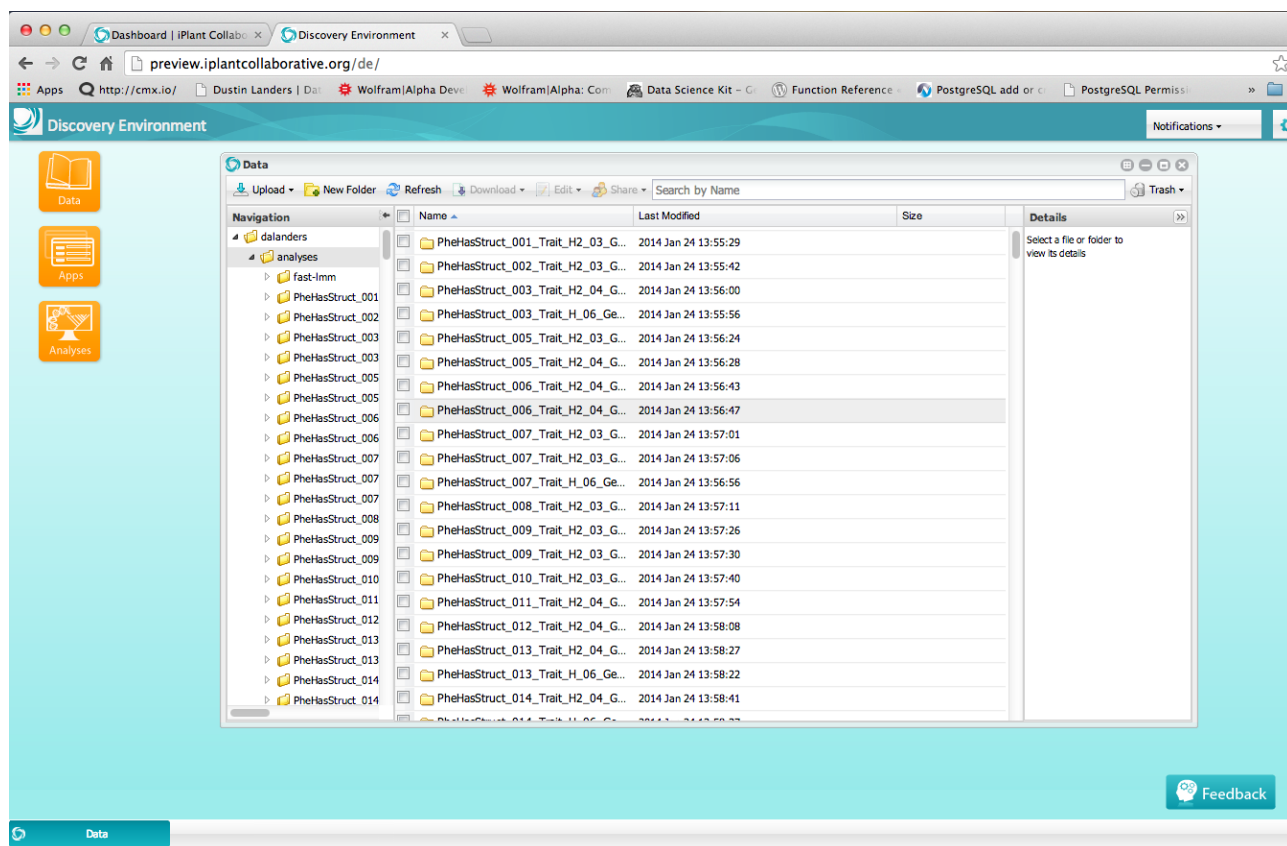

→ If you are testing for effect sizes that are different for different levels of heritability, you may want to create more folders in this step. Actually, for the rest of this process, you would have to everything multiple times for each folder you create here. Say for example, I'm testing three different heritability levels, which will have three different groups of actual effects. Then I would need to run Validate three different times, and thus need to make three folders in this step and separate them that way.

→ Note that in rPlant, you don't have to type your username and you begin iPlant data store file locations immediately. For example, in rPlant, the same location would just be written as "analyses"

Step 2-B) In the DE, create a new folder for your Validate analysis.

The *.qassoc* files in this case, are the files containing the actual analysis that we will use for evaluating the effectiveness of PLINK. Particularly, those are the predicted values for each SNP based on the PLINK algorithm. Now, it's time to compare those predicted values with what we know to be the actual values.

You need to create a folder that we will use for our Validate analysis. In the DE, create a folder called *my_Validate_analysis*.

The goal is to look inside each of the 300 folders we created, and pull out the *.qassoc* file and move it to the folder marked *my_Validate_analysis*.


Since this would take hours, we suggest you launch the application we created called *Aggregate* that provides a fairly easy way to do this (at least, comparatively easier).

Step 2-C) Launch Aggregate.

Once you launched *Aggregate* you will see there is a place to enter your username and password. Go ahead and enter those.


Step 2-D) Enter username and password.

Now, you need to type the out the full location of the folder starting with username where your analyses are located. In this case, all of our PLINK analyses folders are located in the folder *dalanders/analyses*.

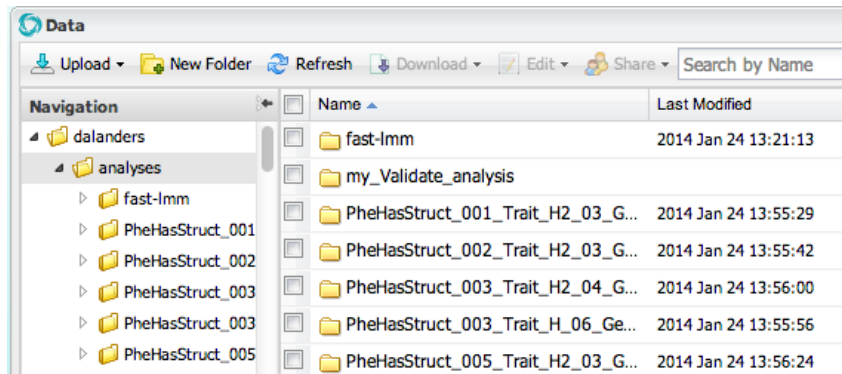

Step 2-E) Enter the folder location where all your tool analyses are located and click *View Folders in Directory*.

We will enter *dalanders/analyses* and click the button *View Folders in Directory*. This will provide us with a list of all the files and folders in that directory on the iPlant data store in the left-most listbox.

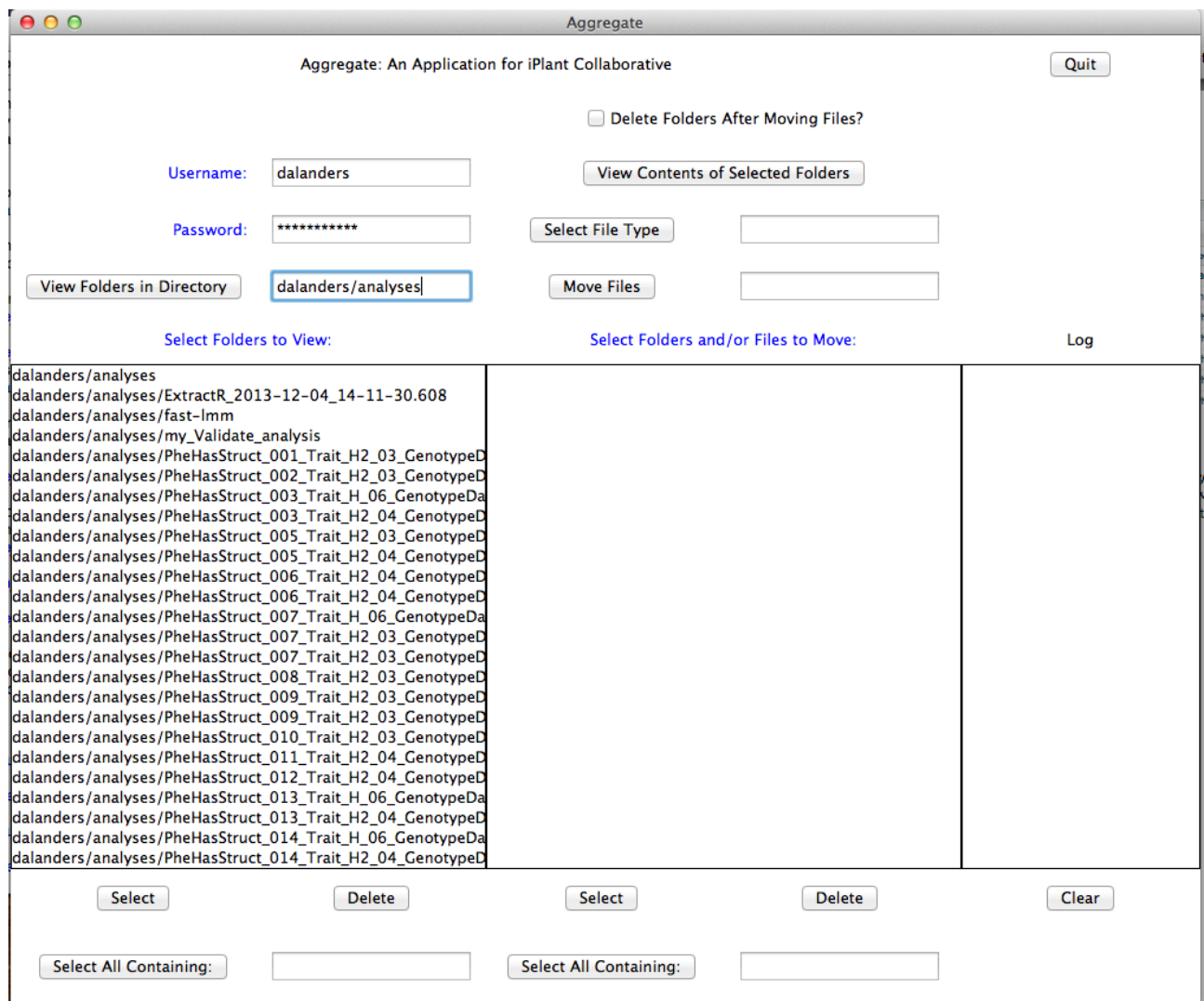

Step 2-F) Select only folders from the known-truth data sets runs.

Now, we want to just select the folders that we just created. You can select them all using Control + Click to select individual ones, or if you wish to just select the whole list you can select all between two files by using Shift + Click.

We also provide a method selecting many at us by using the *Select All Containing:* logic button at the bottom of the listbox.

For example, our analyses folders all contain the words "GenotypeData" or even "Trait".

We can type in the box just to the right of *Select All Containing:* "Trait" and then click the button. This should select just the folders we are interesting obtaining the contents of.

We can scroll up and down to make sure that we have only selected the ones that we need. If we have selected a few that we didn't intend to. We can select those using the mouse and then click the *Delete* button to remove those from the list.

Aggregate: An Application for iPlant Collaborative

Quit

☐ Delete Folders After Moving Files?

Username: dalanders View Contents of Selected Folders

Password: ***** Select File Type

View Folders in Directory dalanders/analyses Move Files

Select Folders to View: Select Folders and/or Files to Move: Log

dalanders/analyses/PheNPStruct_086_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_087_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_088_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_089_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_089_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_090_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_090_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_091_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_091_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_091_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_092_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_092_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_094_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_094_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_094_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_096_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_096_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_096_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_097_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_097_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_098_Trait_H_06_GenotypeDat
dalanders/analyses/PheNPStruct_098_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_098_Trait_H2_04_GenotypeDa
dalanders/analyses/PheNPStruct_099_Trait_H2_03_GenotypeDa
dalanders/analyses/PheNPStruct_100_Trait_H2_03_GenotypeDa

Select Delete Select Delete Clear

Select All Containing: Trait Select All Containing:

Step 2-G) Click *View Contents of Selected Folders*

Once you've got all your folders in the first listbox, click the button on the

right call *View Contents of Selected Folders* and this will iterate over the folders in the left-most listbox and add the contents of all those folders to the right-most listbox. It should take a few minutes, so give it some time. Once its done, see if you can spot the files that you were interested in moving. In our case, they are the files ending in *.qassoc*.

As you can see, every single item in the right-most listbox is a file or folder contained with *one* of the folders in the left-most listbox. Now, we need to select all the files from the right-most listbox that we wish to move to the folder we created in order to run *Validate*.

Step 2-H) Select just the tool outputs from the right-most listbox.

In order to select just the right files that you intend to move, you can use the same *Select All Containing:* logic button we used earlier, or if they all have the same file extension (as ours do with *.qassoc*) then we can use another option.

Above the right-most listbox is a button labelled *Select File Type*. Next to that, we type *qassoc* in and then click the button in order to just select the file types that we want to move.

Aggregate
Aggregate: An Application for iPlant Collaborative
Quit
☐ Delete Folders After Moving Files?
Username: dalanders View Contents of Selected Folders
Password: ***** Select File Type qassoc
View Folders in Directory dalanders/analyses Move Files
Select Folders to View: Select Folders and/or Files to Move: Log

dalanders/analyses/PheNPStruct_086_Trait_H2_04_GenotypeData_001_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_087_Trait_H2_03_GenotypeData_002_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_088_Trait_H2_04_GenotypeData_003_Trait_H_06_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_089_Trait_H2_03_GenotypeData_003_Trait_H2_04_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_089_Trait_H2_03_GenotypeData_005_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_090_Trait_H_06_GenotypeData_005_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_090_Trait_H2_03_GenotypeData_006_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_091_Trait_H_06_GenotypeData_006_Trait_H2_04_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_091_Trait_H2_03_GenotypeData_007_Trait_H_06_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_091_Trait_H2_04_GenotypeData_007_Trait_H2_03_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_092_Trait_H_06_GenotypeData_007_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_092_Trait_H2_04_GenotypeData_008_Trait_H2_03_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_094_Trait_H_06_GenotypeData_009_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_094_Trait_H2_03_GenotypeData_009_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_094_Trait_H2_04_GenotypeData_010_Trait_H_06_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_096_Trait_H_06_GenotypeData_011_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_096_Trait_H_06_GenotypeData_012_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_096_Trait_H2_04_GenotypeData_013_Trait_H_06_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_097_Trait_H_06_GenotypeData_013_Trait_H2_04_GenotypeData_NoMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_097_Trait_H2_03_GenotypeData_014_Trait_H_06_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_098_Trait_H_06_GenotypeData_014_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_098_Trait_H2_03_GenotypeData_015_Trait_H_06_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_098_Trait_H2_04_GenotypeData_015_Trait_H2_04_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_099_Trait_H2_03_GenotypeData_016_Trait_H2_04_GenotypeData_8pctMissing_--assoc.qassoc	View:success
dalanders/analyses/PheNPStruct_100_Trait_H2_03_GenotypeData_017_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	View:success

Select Delete Select Delete Clear

Select All Containing: Trait Select All Containing:

Step 2-I) Move the files to your *Validate* analysis folder.

Now, we can move the files. In the text box next to the *Move Files* button enter the full path of the folder we created just for this purpose in Step 2-B. In our case, we named it *my_Validate.analysis* so we would enter in to the text box: *dalanders/analyses/my_Validate.analysis*. Then click *Move Files* in order to begin iterating over that list and sending the *.qassoc* files to the folder you created.

This process can take about as long or longer than it did to actually view the files. Keep in mind, that this process uses the iPlant Foundation API system, so for each item in the right-most list, a request is send to the API to move that file to the folder you specified.

It may also be worth noting that by the time you read this manual, this could be a fully-integrated step in the iPlant Discovery Environment, as there was talk about this functionality being available given more time. However, we found it necessary, even though this part of the problem is almost entirely logistical, to be able to select files from multiple files and move them in to a single aggregate folder.

Aggregate

Aggregate: An Application for iPlant Collaborative

Quit

☐ Delete Folders After Moving Files?

Username: dalanders

View Contents of Selected Folders

Password: *****

Select File Type

qassoc

View Folders in Directory

dalanders/analyses

Move Files

ses/my Validate analysis

Select Folders to View:

Select Folders and/or Files to Move:

Log

dalanders/analyses/PheHasStruct_019_Trait_H_06_GenotypeData_001_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_019_Trait_H2_04_GenotypeData_002_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_020_Trait_H2_03_GenotypeData_003_Trait_H_06_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_021_Trait_H_06_GenotypeData_003_Trait_H2_04_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_021_Trait_H2_03_GenotypeData_005_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_022_Trait_H2_03_GenotypeData_005_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_022_Trait_H2_04_GenotypeData_006_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_023_Trait_H2_03_GenotypeData_006_Trait_H2_04_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_024_Trait_H_06_GenotypeData_07_Trait_H_06_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_025_Trait_H2_03_GenotypeData_007_Trait_H2_03_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_025_Trait_H2_03_GenotypeData_007_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_025_Trait_H2_04_GenotypeData_008_Trait_H2_03_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_026_Trait_H2_03_GenotypeData_009_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_027_Trait_H2_04_GenotypeData_009_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_028_Trait_H_06_GenotypeData_010_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_028_Trait_H2_03_GenotypeData_011_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_029_Trait_H_06_GenotypeData_012_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_030_Trait_H_06_GenotypeData_013_Trait_H_06_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_030_Trait_H2_04_GenotypeData_013_Trait_H2_04_GenotypeData_NoMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_032_Trait_H_06_GenotypeData_014_Trait_H_06_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_032_Trait_H_06_GenotypeData_014_Trait_H2_04_GenotypeData_3pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_033_Trait_H_06_GenotypeData_015_Trait_H_06_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_033_Trait_H2_03_GenotypeData_015_Trait_H2_04_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_034_Trait_H_06_GenotypeData_016_Trait_H2_04_GenotypeData_8pctMissing_--assoc.qassoc	
dalanders/analyses/PheHasStruct_036_Trait_H2_03_GenotypeData_017_Trait_H2_03_GenotypeData_8pctMissing_--assoc.qassoc	

Select

Delete

Select

Delete

Clear

Select All Containing:

Trait

Select All Containing:

3.2 How to Use the Aggregate Function in rPlant

3.3 How to Use Aggregate on Atmosphere

4 Validating with Validate

4.1 How to Use Validate on the DE

Step 3-A) First, take a look at one of your analyses documents on the DE. Check for four things: The nature of the column delimitation, the name of the SNP column, the name of P-value or other importance scoring column, and the name of an effect size column if there is one.

For example, our SNP column is SNP, our importance column is P, and our effect size column is BETA. You are looking for those four things mentioned.


The goal is to assess how well PLINK was at identifying those known SNPs as compared to the ones where there are no effects. We are also going to include a SNP effects text file, so that we can measure the effectiveness for PLINK at determining those.

The screenshot shows a web-based data interface. On the left is a file explorer with a navigation pane showing a directory structure: 'dalanders' > 'analyses' > 'fast-Imm' > 'my_Val' > 'PheHas'. The main pane shows a list of files: 'rait_H2_03_GenotypeData_NoMissing.map', 'rait_H2_03_GenotypeData_NoMissing.ped', and 'PheHasStruct_001_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc'. The file 'PheHasStruct_001_Trait_H2_03_GenotypeData_NoMissing_--assoc.qassoc' is selected, and its details are shown on the right: 'Last Modified : 2014 Jan 24', 'Date Submitted : 2014 Jan 24', and 'Permissions : Owner'. Below the file list is a table with 9 columns: CHR, SNP, BP, NMISS, BETA, SE, R2, T, and P. The table contains 24 rows of data for SNPs on chromosome 1, labeled AWA1 through AWA24. The data includes values for BETA, SE, R2, T, and P for each SNP.


CHR	SNP	BP	NMISS	BETA	SE	R2	T	P
1	AWA1	1	512	0.5428	0.321	0.005573	1.691	0.09151
1	AWA2	2	512	1.373	2.378	0.0006536	0.5776	0.5638
1	AWA3	3	512	-3.738	2.373	0.004842	-1.575	0.1158
1	AWA4	4	512	-3.738	2.373	0.004842	-1.575	0.1158
1	AWA5	5	512	2.264	0.9708	0.01056	2.333	0.02006
1	AWA6	6	512	0.9807	1.067	0.001653	0.919	0.3585
1	AWA7	7	512	3.003	1.678	0.006239	1.789	0.07415
1	AWA8	8	512	1.607	0.796	0.007929	2.019	0.04401
1	AWA9	9	512	1.411	1.191	0.002744	1.185	0.2367
1	AWA10	10	512	0.7191	1.068	0.0008889	0.6736	0.5009
1	AWA11	11	512	1.04	0.6927	0.004404	1.502	0.1337
1	AWA12	12	512	1.32	1.066	0.002997	1.238	0.2162
1	AWA13	13	512	0.9997	0.3982	0.01221	2.511	0.01236
1	AWA14	14	512	-2.853	2.375	0.00282	-1.201	0.2303
1	AWA15	15	512	0.2873	0.9043	0.0001978	0.3176	0.7509
1	AWA16	16	512	3.042	1.185	0.01275	2.566	0.01056
1	AWA17	17	512	1.809	0.9726	0.006735	1.86	0.06352
1	AWA18	18	512	1.704	0.973	0.005978	1.751	0.08048
1	AWA19	19	512	0.9714	0.6426	0.004461	1.512	0.1312
1	AWA20	20	512	-2.045	2.377	0.001449	-0.8604	0.39
1	AWA21	21	512	0.04449	0.3365	3.427e-05	0.1322	0.8949
1	AWA22	22	512	-2.245	2.377	0.001746	-0.9445	0.3454
1	AWA23	23	512	0.494	0.7989	0.000749	0.6183	0.5367
1	AWA24	24	512	-0.1613	2.379	9.02e-06	-0.06783	0.946

Step 3-B) Gather known SNP and known effect text files.

The known SNP file contains essentially just text with a list of all the SNPs with known effects in our case, its a list of 35 SNPs. For us, these are called syntruth.txt and synbetas03.txt respectively and we have included screenshots so you can get an idea of what they need to look like.

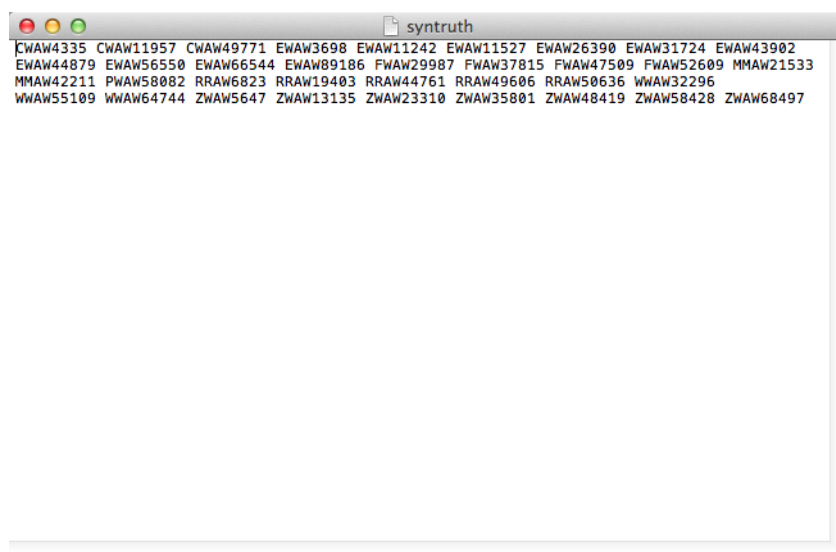

The known effects text file is the same way and it *must* be in the same order as the SNPs. This means (looking at the above screenshots) that SNP CWAW1335 has effect 0.0324 and so on. This information ordering information is used by *Validate* to assign true effect sizes to all the SNPs in your outputs (and in doing so labels all others as zero). Also, *Validate* uses the information from the known SNPs file to label all SNPs in your output as either TRUE or FALSE based on whether or not they are included in that file.

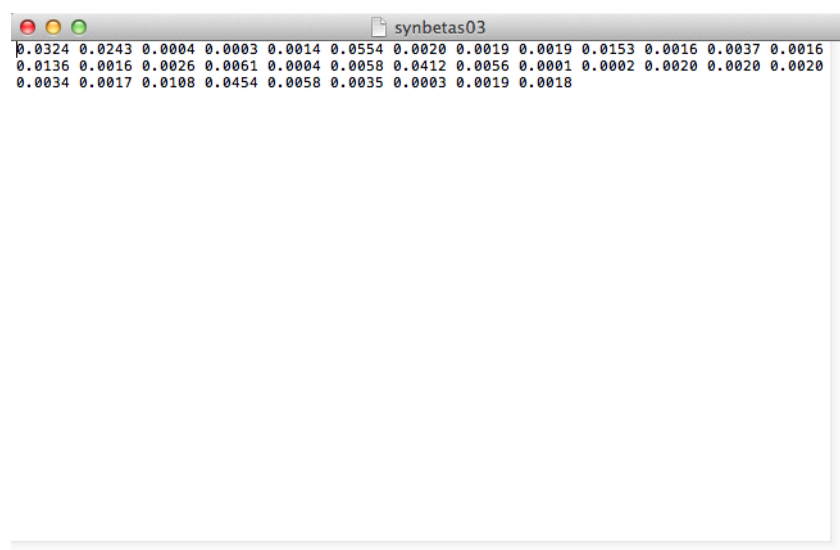

The screenshot shows a text file named 'synbetas03' with a grid of numerical values. The values are arranged in 3 rows and 14 columns. The first row contains values: 0.0324, 0.0243, 0.0004, 0.0003, 0.0014, 0.0554, 0.0020, 0.0019, 0.0019, 0.0153, 0.0016, 0.0037, 0.0016, 0.0016. The second row contains values: 0.0136, 0.0016, 0.0026, 0.0061, 0.0004, 0.0058, 0.0412, 0.0056, 0.0001, 0.0002, 0.0020, 0.0020, 0.0020, 0.0020. The third row contains values: 0.0034, 0.0017, 0.0108, 0.0454, 0.0058, 0.0035, 0.0003, 0.0019, 0.0018, and the rest of the cells are empty.

This is how performance measures are inherently generated, so this *must* be included. By default, *Validate* includes classification measures such as AUC and so you must include known SNPs file. However, the known effects file is optional.

Step 3-C) Upload your known SNPs file and the known effects file using the DE.

Upload both of these text files using the Discovery Environment. This is extremely easy to do. It also does not matter where you upload them to.

Step 3-D) Run Validate!

It's time to run *Validate* finally! Open up the Apps icon on the DE and select *Validate*. Then you need to just fill in the form. Select a name and a description of your choice.

- 5 Visualizing the Validation**
 - 5.1 How to Use the Demonstrate R Package**
 - 5.2 How to Use Demonstrate on the DE**
 - 5.3 How to Use Demonstrate on Atmosphere**
- 6 Acknowledgements**
- 7 How to Interpret Performance Measures**
 - 7.1 AUC**
 - 7.2 H**
 - 7.3 KS**
 - 7.4 TPR**
 - 7.5 FPR**
 - 7.6 Accuracy**
 - 7.7 MAE**
 - 7.8 RMSE**

hmeasure.net